



HEK293-BFP Stable Cell Line Manual

Catalog Number	Product name	Amount
SC102	HEK293- BFP (Bsd) stable cells	1 ml / vial in 80% DMEM, 10% FBS, 10% DMSO


Storage:

Upon received, place vial in Liquid Nitrogen for long-term storage, or saved in -80oC for short-time storage up to 1 week.


1. Product Description

The HEK293 Cell Line is a permanent line established from primary embryonal human kidney cells transformed with sheared human adenovirus type 5 DNA. The expressed E1A adenovirus gene allows these cells to produce very high levels of protein.

HEK293-BFP cells were transformed from the HEK293 cell line and stably express the human-codon, evergreen Blue Fluorescent Protein (**BFP**), wavelength peak: 365-380nm (Ex) / 450-470nm (Em). The cell line was established by transduction with BFP expression lentivirus containing a Blasticidin-resistance gene for selection and maintenance of the stable cells. BFP is constitutively expressed at high-levels under the CMV promoter. The following expression construct was integrated into cell's genome. Each cell demonstrates strong fluorescent signal under microscope (see image below).

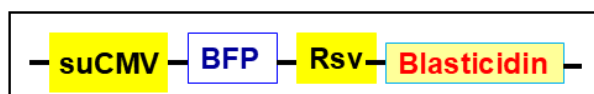

CAT#: **SC102**

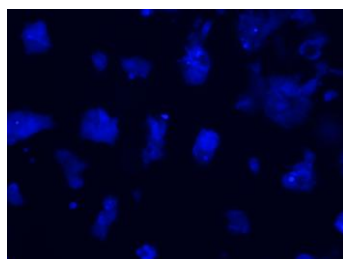

BFP Filter



Bright Field



2. Culture procedures

1. Thaw the frozen vial of cells quickly in a 37°C water bath (1~3min), decontaminate the outside of the vial with 70% ethanol.
2. Transfer the entire contents of the cryovial into a T-75 cm² flask containing 20 ml of pre-warmed complete medium. Incubate the cells overnight in a 37°C incubator, 5% CO₂.
3. On the following day, replace the medium with 20 ml of prewarmed, complete medium.
4. Incubate the cells and monitor cell density.
5. Pass cells (1:5 to 1:10 dilution) using 0.25% Trypsin-EDTA solution when the culture reaches ~90% confluent.
6. Freeze cells at a density of $\sim 3 \times 10^6$ cells/ml using 90% complete medium with 10% DMSO.

3. Complete medium

DMEM (high glucose)
2mM L-glutamine
10% Fetal Bovine Serum (FBS)
0.1 mM MEM Non-Essential Amino Acids (NEAA)
1% Pen-strep or Antibiotic-antimycoplasma

4. Quality Control

Each vial contains $\sim 2 \times 10^6$ cells with >95% viability before freezing. Cells are verified to be free of bacteria, viruses, and mycoplasma.

5. Warranty and user terms

- 1) This product is warranted to perform as described when used in accordance with this manual. GenTarget's sole remedy for breach of warranty should be, at the option of GenTarget, to repair or replace the product if this product does not meet the stated quality standard.
- 2) By paying the purchase price, the buyer is granted a non-transferable, non-exclusive license to use the product. This product is sold **for research and development purposes only**. For commercial use, please contact Gentarget Inc for license.
- 3) This product is limited to the laboratory that the product is delivered to. This Product is not for resale, distribution, or transfer for any purpose, including transfer of the Product as a component of any product(s);



GenTarget Inc

7930 Arjons Drive, Suite B
San Diego, CA 92126, USA
Phone: 1(858) 2656446
Email: orders@gentarget.com
Fax: 1(800) 3804198

GenTarget will retain all rights for this Product's license and other intellectual property.

- 4) This Product should be used only for non-profit purposes including any products and services usages; furthermore, **research use only** means that this product is excluded, without limitation, from resale, repackaging, or modification for the making or selling of any commercial product(s) or service(s) without the written approval of GenTarget. You may contact our Business Development department at support@gentarget.com for product proprietary information.
- 5) GenTarget is not liable, and does not have any responsibility or liability, whatsoever for any direct and indirect, consequential, or other damages resulting from using this Product.
- 6) Gentarget **do not** provide the protected reporter's sequences information for all our cell line products unless a license is purchased.